

SINGLE USE



S4004

Our Four Piece Restorative Kit includes our Ball Burnisher, Flat Plastic, London College Tweezers and with the inclusion of our small spoon excavator, the S4004 is perfect for paediatric restorations.

Single Use Dental Instruments Pty Ltd

1/24 Hoopers Road, Kunda Park
QLD 4556 Australia
sales@sudinstruments.com
07 5445 6486

Please see website link for kit details [here](#).

ERGONOMIC • PRECISE • QUALITY • STERILE • SINGLE USE

BB56

Ball Burnisher


Description:

The perfect blending tool. Our ball burnisher has both small and large tip ends to effectively condense and shape the composite against the cavity wall surfaces.

Features:

- Double-ended instrument
- Precision manufactured stainless steel tips
- Ergonomic handle with ridges and grooves for improved grip
- Stainless steel stem for added stability
- Latex-free
- Arrives sterile and ready for use
- CS1 Data matrix
- CE and TGA Certified

Materials: Stainless steel, Acrylonitrile butadiene styrene (ABS)

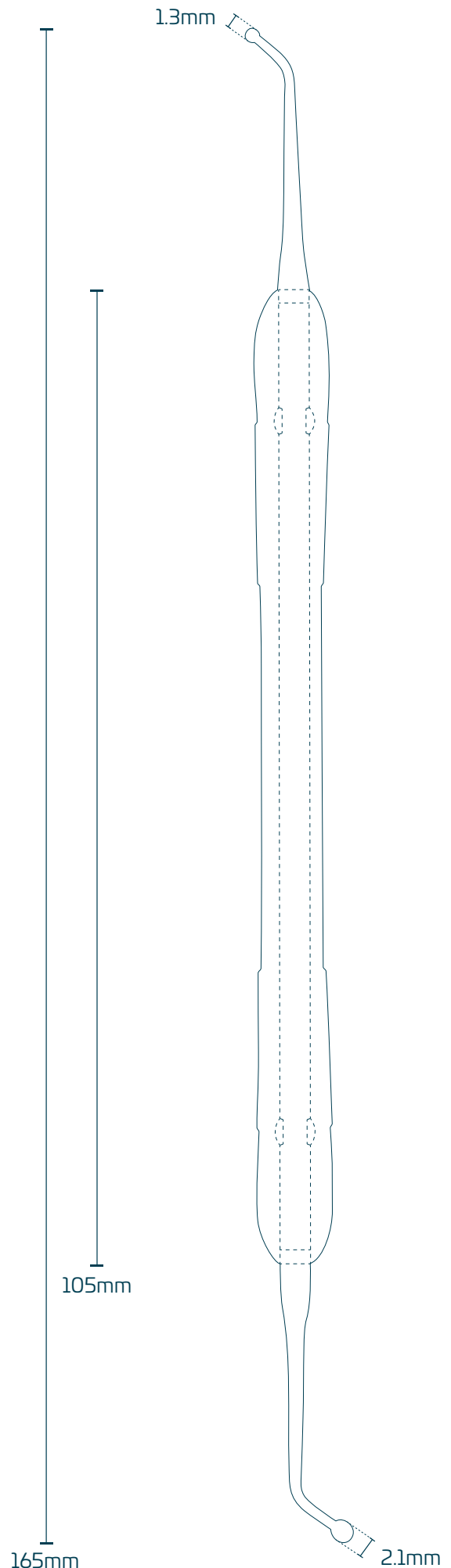

XC67

Small Spoon Excavator (1.2mm)


Description:

One of our bestsellers, our Small Spoon Excavator is a popular choice for removing caries within deep cavities.

Features:

- Double-ended instrument
- Precision manufactured stainless steel tips
- Small spoon diameter; 1.2mm
- Ergonomic handle with ridges and grooves for improved grip
- Stainless steel stem for added stability
- Latex-free
- Arrives sterile and ready for use
- CS1 Data matrix
- CE and TGA Certified

Materials: Stainless steel, Acrylonitrile butadiene styrene (ABS)

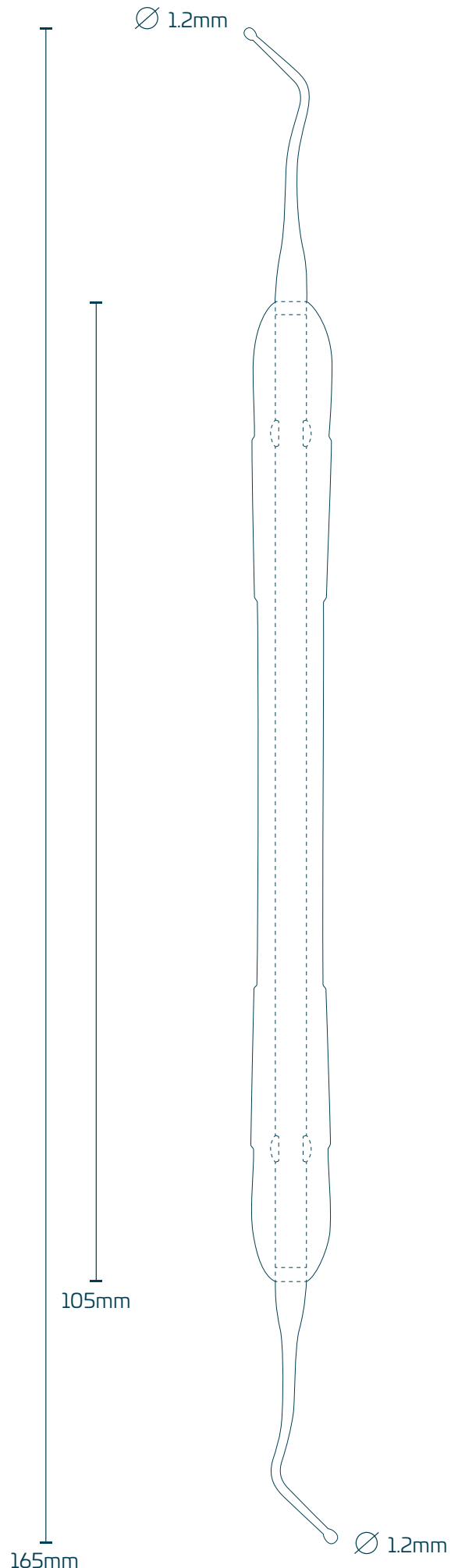

FP54

Flat Plastic


Description:

Our double ended Flat Plastic features long, flat tips that are used to carry and shape composite materials.

Features:

- Double-ended instrument
- Precision manufactured stainless steel tips
- Tip diameters; 2.2mm/1.8mm
- Ergonomic handle with ridges and grooves for improved grip
- Stainless steel stem for added stability
- Latex-free
- Arrives sterile and ready for use
- CS1 Data matrix
- CE and TGA Certified

Materials: Stainless steel, Acrylonitrile butadiene styrene (ABS)

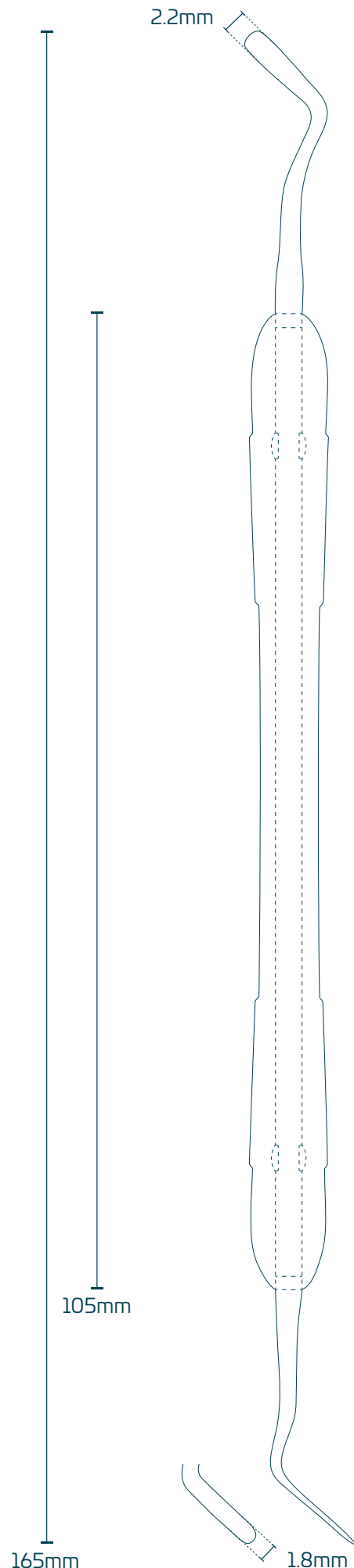

TW82

London College Tweezers


Description:

Strong and sturdy with enhanced grip features on the handle and tips. Our premium London College Tweezers are an essential during procedures.

Features:

- Ergonomic handle
- Serrated non-slip tips
- Latex-free
- Arrives sterile and ready for use
- GS1 Data matrix
- CE and TGA Certified

Materials: Stainless steel

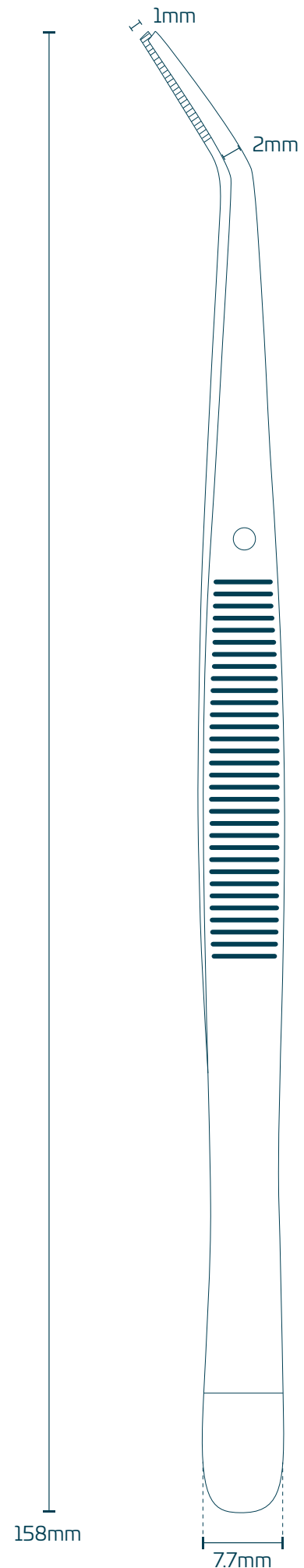

Our Products

Our unique, single-use instruments were designed with simplicity in mind. Featuring a broad variety of diagnostic and restorative instruments, our goal is to remove daily stress for clinicians and improve patient safety. Not only convenient, we have created a product range that is focused on quality and sustainability. Our single use instruments arrive sterile and ready for immediate use. The ergonomic handle is made from ABS plastic and features ridges and grooves for an improved grip. The stainless-steel stem provides additional weight and durability to ensure it is fit for purpose. Precision manufactured stainless steel tips replicate the design and function of re-usable instruments.

Packaging and materials

We are committed to reducing the environmental impact of our products. As such our products and packaging are manufactured with recyclable materials. Additionally, we source recycled cardboard to create our benchtop dispensers and shipping containers.

Quality assurance

All our instrument manufacturing has CE, TGA and ISO certification to ensure quality compliance with industry standards.

Tyvek® medical dialysing strip

A Tyvek® medical dialysing strip is incorporated into all our product sleeves. Comprising of flashspun and bonded high-density polyethylene (HDPE) filaments, the strip provides a microbial barrier and exceptional tear and puncture resistance to maintain product integrity during transit.

