

# SINGLE USE



**This kit also features:**

2 Face Masks with Ear Loops  
Patient Bib with Ties  
2x Cotton Roll No 2  
Tissue Wipe  
Tray

## D3550

Convenient for mobile dentists and domiciliary care. With our Front Face Mirror, Perio Probe/Sickle and Examination Tweezers at your fingertips, our diagnostic pack is making patient examinations quick and easy. In addition, we also have 2x Face Masks with Ear Loops, 2x Cotton Rolls No.2, Patient Bib with Ties, Tissue wipe & Tray for your convenience.

Single Use Dental Instruments Pty Ltd

1/24 Hoopers Road, Kunda Park  
QLD 4556 Australia  
sales@sudinstruments.com  
07 5445 6486

ERGONOMIC • PRECISE • QUALITY • STERILE • SINGLE USE

# FF53

## Front Face Mirror

### Description:


An essential item in any dental practice our Front Face Mirror is new and improved. With a focus on comfort, our ergonomic handle features additional ridges and grooves to improve grip stability. With the addition of a zinc alloy tip for percussion testing and an anti-fog mirror surface.

### Features:

- Zinc alloy tip for percussion testing
- Ergonomic handle with ridges and grooves for improved grip
- Latex-free
- Arrives sterile and ready for use
- GS1 Data Matrix
- CE and TGA Certified

Materials: Zinc Alloy, Acrylonitrile butadiene styrene (ABS), Polycarbonate (PC), Aluminiumbutadiene styrene (ABS)

**SINGLE USE**

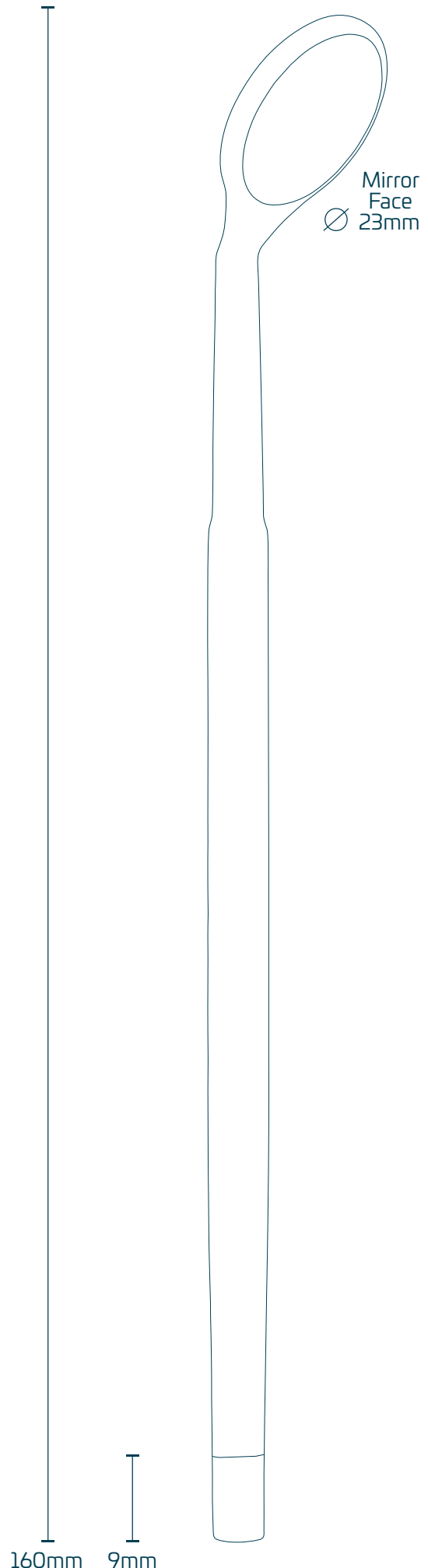

# PB84

## No. 8/CPITN-C Probe

### Description:


Ensuring that our product range caters for everyone, our PB84 is popular amongst dentists trained within the UK.

### Features:

- Combination instrument
- Precision manufactured stainless steel tips
- Ergonomic handle with ridges and grooves for improved grip
- Calibrations; CPITN-C: 3.5/5.5/8.5/11.5mm plus 0.5mm wide ball tip
- Latex-free
- Arrives sterile and ready for use
- GS1 Data matrix
- CE and TGA Certified

Materials: Stainless steel, Acrylonitrile butadiene styrene (ABS)

**SINGLE USE**

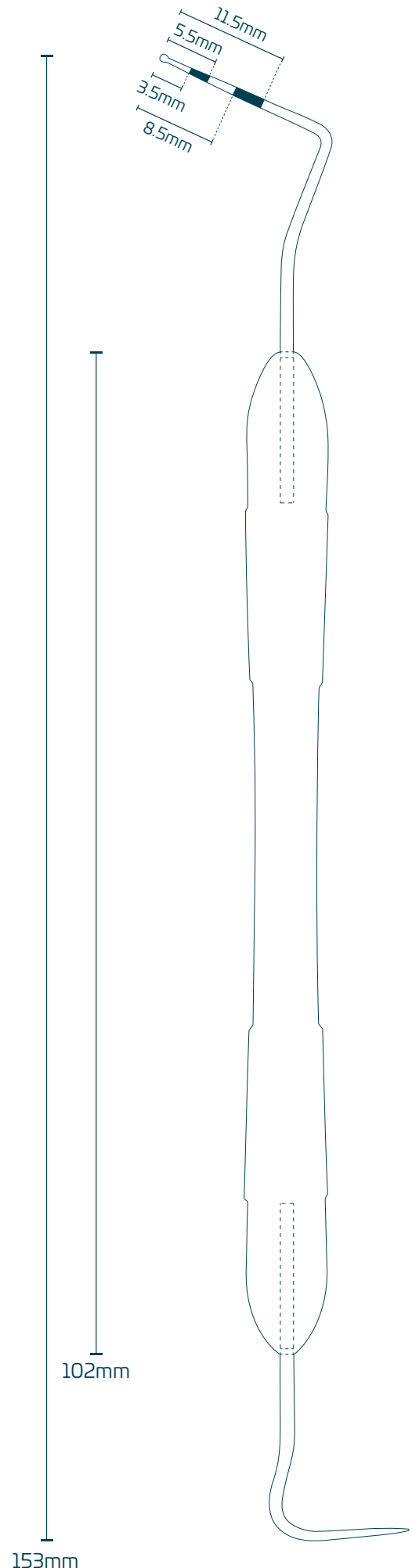

# TW81

## Examination Tweezers


### Description:

Our economy Tweezer is a popular choice for both diagnostic and restorative consultations. Featuring serrated, non-slip tips to ensure you have complete control during procedures.

### Features:

- Ergonomic handle
- Serrated non-slip tips
- Latex-free
- Arrives sterile and ready for use
- GS1 Data matrix
- CE and TGA Certified

Materials: Steel plate, zinc coating

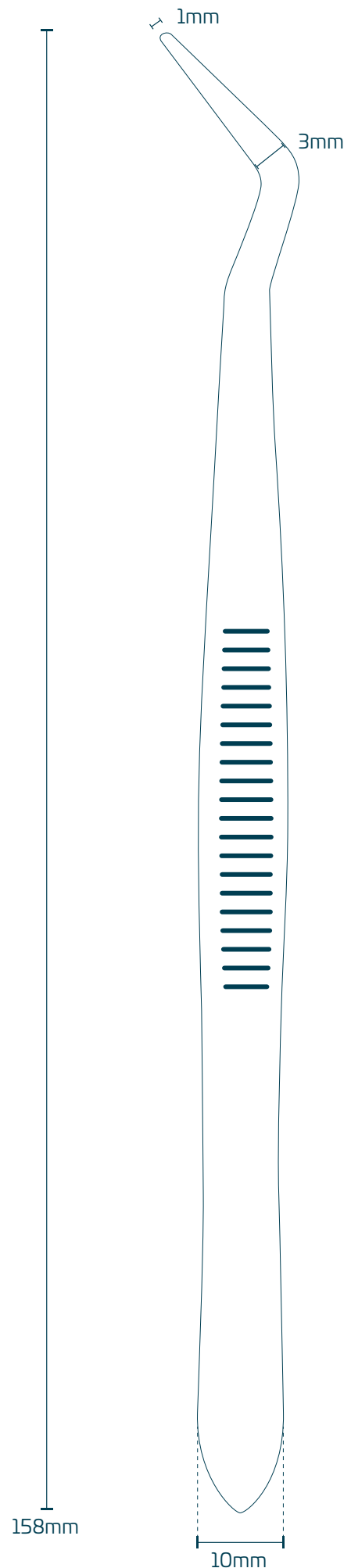

**SINGLE USE**

# Our Products

Our unique, single-use instruments were designed with simplicity in mind. Featuring a broad variety of diagnostic and restorative instruments, our goal is to remove daily stress for clinicians and improve patient safety. Not only convenient, we have created a product range that is focused on quality and sustainability. Our single use instruments arrive sterile and ready for immediate use. The ergonomic handle is made from ABS plastic and features ridges and grooves for an improved grip. The stainless-steel stem provides additional weight and durability to ensure it is fit for purpose. Precision manufactured stainless steel tips replicate the design and function of re-usable instruments.

## Packaging and materials

We are committed to reducing the environmental impact of our products. As such our products and packaging are manufactured with recyclable materials. Additionally, we source recycled cardboard to create our benchtop dispensers and shipping containers.

## Quality assurance

All our instrument manufacturing has CE, TGA and ISO certification to ensure quality compliance with industry standards.

## Tyvek® medical dialysing strip

A Tyvek® medical dialysing strip is incorporated into all our product sleeves. Comprising of flashspun and bonded high-density polyethylene (HDPE) filaments, the strip provides a microbial barrier and exceptional tear and puncture resistance to maintain product integrity during transit.

