


SINGLE USE



D3003

**Three Piece |
Orthodontic Kit**

FF53 | Front Face Mirror with Percussion Tip

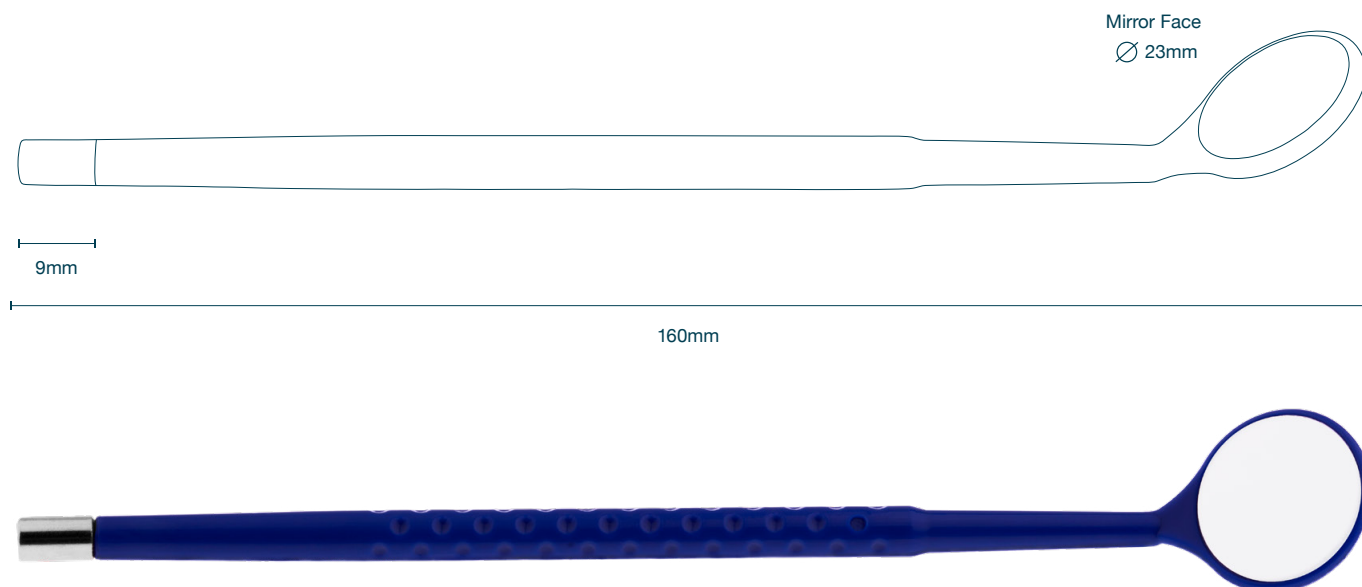

Description

An essential item in any dental practice our Front Face Mirror is new and improved. With a focus on comfort, our ergonomic handle features additional ridges and grooves to improve grip stability. With the addition of a zinc alloy tip for percussion testing and an anti-fog mirror surface.


Features

- Zinc alloy tip for percussion testing
- Ergonomic handle with ridges and grooves for improved grip
- Latex-free
- Arrives sterile and ready for use
- GS1 Data Matrix
- CE Certified

Materials

- Zinc Alloy
- Acrylonitrile butadiene styrene (ABS)
- Polycarbonate (PC)
- Aluminiumbutadiene styrene (ABS)

PB89 | Sickle Explorer & Perio Probe

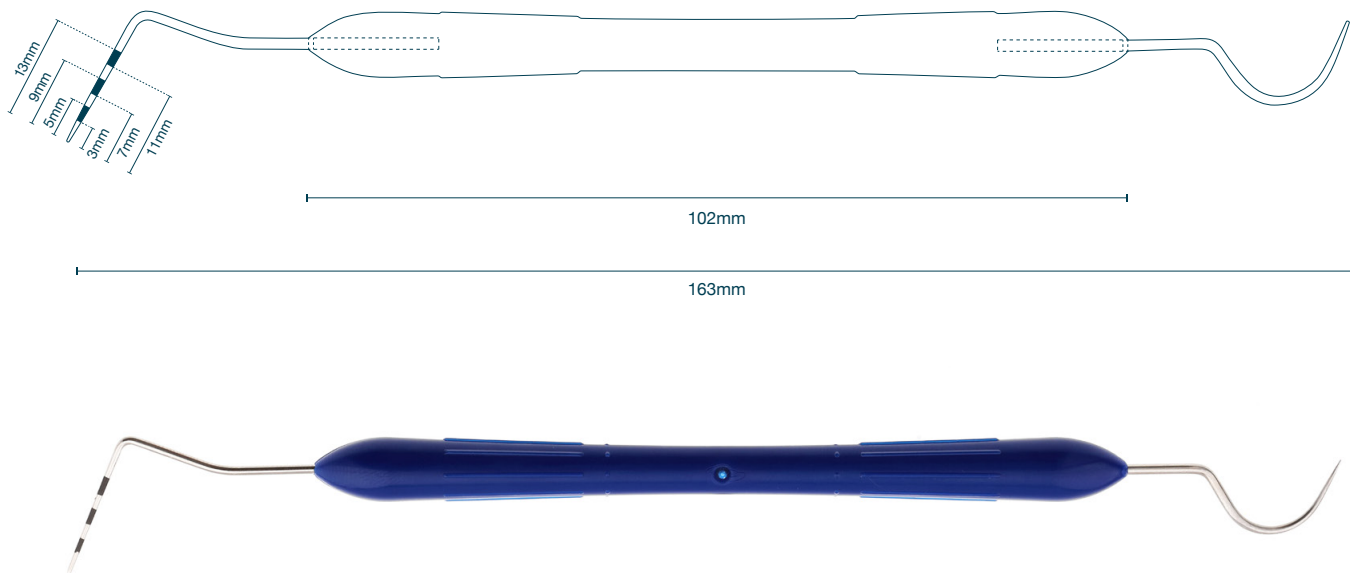

Description

A customer favourite, our perio probe/sickle has easily identifiable markings to assist in delivering accurate assessments. This combination instrument can be purchased individually or within our top selling diagnostic kits; D3001 and D3006.


Features

- Combination instrument
- Precision manufactured stainless steel tips
- Ergonomic handle with ridges and grooves for improved grip
- Calibrations; Perio Probe: 3/5/7/9/11/13mm
- Latex-free
- Arrives sterile and ready for use
- GS1 Data matrix
- CE Certified

Materials

- Stainless steel
- Acrylonitrile butadiene styrene (ABS)

Orthodontic Ruler

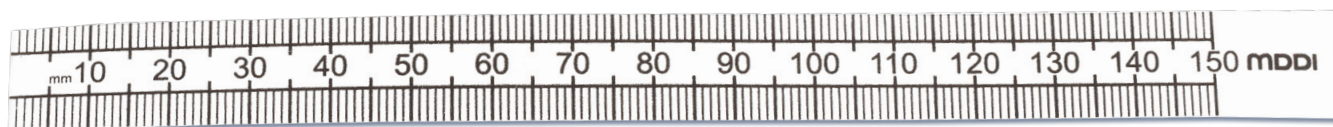

Description

Sterile, orthodontic ruler measuring 150mm (1mm increments). Flexible polypropylene material allowing for accurate intra-oral measurements.

Features


- Flexible Design, Bends without breaking for easy maneuverability.
- Zero-Start Scale, Measurements begin right at the edge for precision accuracy.

Materials

- Polypropylene (with ink print)

Products arrive in a protective, sterile sleeve allowing for immediate use.


Removable label with GS1 Data Matrix and product code for electronic data recording and patient logbooks.

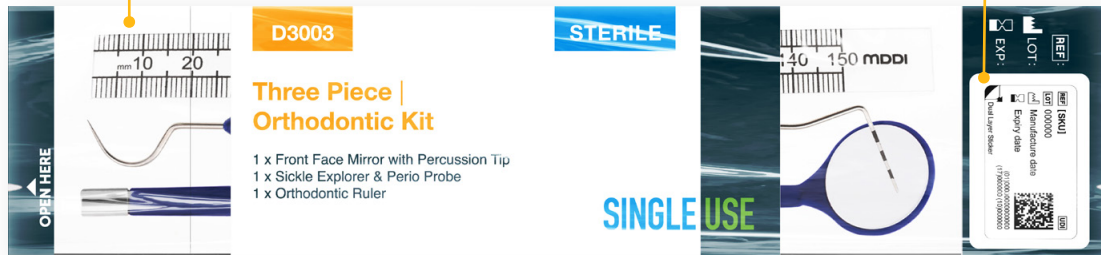

Easy access tear design for aseptic instrument opening.

Identifiable instructions for use and product storage requirements.

Sterilised with Ethylene Oxide.

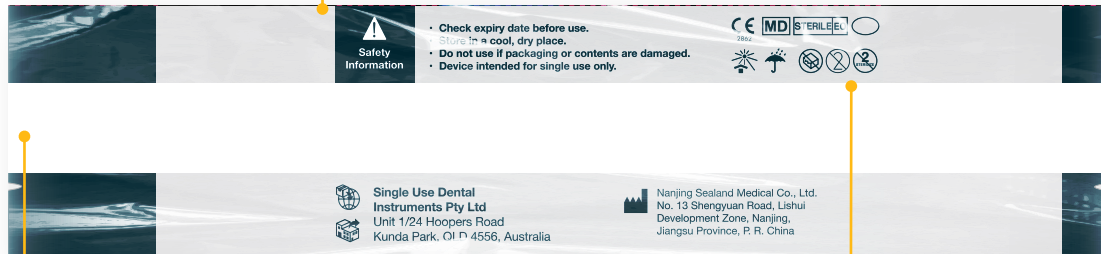

Tyvek® medical dialysing strip.

CE & regulatory identification marks in line with industry standards.

Our Products

Our unique, single-use instruments were designed with simplicity in mind. Featuring a broad variety of diagnostic and restorative instruments, our goal is to remove daily stress for clinicians and improve patient safety. Not only convenient, we have created a product range that is focused on quality and sustainability. Our single use instruments arrive sterile and ready for immediate use. The ergonomic handle is made from ABS plastic and features ridges and grooves for an improved grip. The stainless-steel stem provides additional weight and durability to ensure it is fit for purpose. Precision manufactured stainless steel tips replicate the design and function of re-usable instruments.

Packaging and materials

We are committed to reducing the environmental impact of our products. As such our products and packaging are manufactured with recyclable materials. Additionally, we source recycled cardboard to create our benchtop dispensers and shipping containers.

Quality assurance

All our instrument manufacturing has CE and ISO certification to ensure quality compliance with industry standards.

Tyvek® medical dialysing strip

A Tyvek® medical dialysing strip is incorporated into all our product sleeves. Comprising of flashspun and bonded high-density polyethylene (HDPE) filaments, the strip provides a microbial barrier and exceptional tear and puncture resistance to maintain product integrity during transit.

SINGLE USE

Single Use Dental Instruments Pty Ltd

07 5445 6486 | singleuse.au
sales@sudinstruments.com

1/24 Hoopers Road
Kunda Park, QLD Australia 4556
